


## THE BEST IN BARIATRIC WEIGHING



**Model 6855**  
**Waist-High Stand-On Scale**  
600 lb x .2 lb /  
270 kg x .1 kg



**Model 6868**  
**Bariatric Flip Seat Scale**  
800 lb x .2 lb /  
300 kg x .1 kg

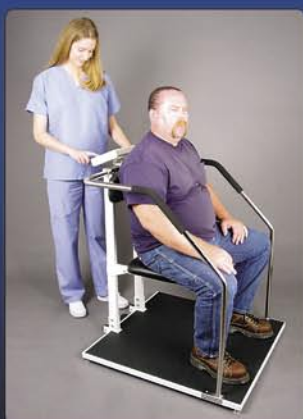

**E**ngineered specifically for weighing obese and unsteady patients, Detecto's 6855 and 6868 portable high-capacity scales offer more features than any other scale in their class. Heavy-duty, yet lightweight, these scales' balanced design makes them easy to move. Both models feature low-profile platforms with non-skid surface and sturdy handrails with padded grips. The unique patent pending design of the 6868 offers an extra-large, 34" x 32" platform with flip seat, enabling patients to be weighed standing or sitting, whichever is most comfortable.

Our digital weight indicator is the industry's most advanced, with easy to use "quick weigh" keys, plus Body Mass Index calculation. The auto-zero function ensures accurate weight each time without having to press the zero key; the display "hold" feature retains weight so you can attend to the patient first, then record the weight. Powered by your choice of either batteries or AC adapter, these scales also have auto-off function to conserve energy.



## Flip-Up

### Seat

*The 6868's unique flip seat adds comfort to the weighing of obese patients and those requiring special assistance. When finished, simply raise it up for those who prefer standing during weighing.*



## Completely Portable

*Both the 6868 and 6855 tip back for easy transporting. The large diameter wheels and handrail make it easy to move from place to place. The compact design permits the scales to fit through doors and into patient's rooms for easy weighing (6868 requires door opening greater than 32").*

Feature	6855	6868
Highly accurate digital weight display with Body Mass Index	•	•
Padded foam black handgrips	•	•
Extra strength wrap-around tubular stainless steel handrail	•	•
LCD weight display can be read from either side of the platform	•	•
Large, heavy duty rubber wheels and tips back for easy transporting	•	•
Skid-proof and soft, but durable black rubber mat	•	•
Completely portable - powered by 6 "C" size Alkaline batteries	•	•
Oversized non-skid platform (34" x 32")		•
Unique flip-up seat allows weighing in a standing or sitting position		•

## Mutual Specifications

Weight Display	.7" (18 mm) high, 7 segment LCD with Body Mass Index
Controls	On/off button, lb/kg Hold/release, Zero Automatic zero
Display Functions	Auto power off Hold/release of display Motion detection
Power	9 VDC - Six "C" size batteries
Color	White with black mat, seat and hand grips
Options	115 VAC, 60 Hz adapter (with 6' power cord for long term stationary operation)

## 6855 Specifications

Capacity	600 lb x .2 lb / 270 kg x .1 kg
Platform Size	18" x 14" / 457 mm x 356 mm
Shipping Weight	70 lb / 32 kg



## 6868 Specifications

Capacity	800 lb x .2 lb / 300 kg x .1 kg
Platform Size	34" x 32" / 864 mm x 813 mm
Shipping Weight	147 lb / 67 kg



Cardinal Scale reserves the right to improve, enhance or modify features and specifications without prior notice.



## Sold By:

**OPS Medical, LLC** An Authorized Detecto Distributor  
**Free: 1(800) 430-0048 Local: 1(410) 553-9293**  
 email [sales@opsmedical.com](mailto:sales@opsmedical.com)  
[www.opsmedical.com](http://www.opsmedical.com)